



Perifusion Chamber

Model No. PERI-CHAMBER

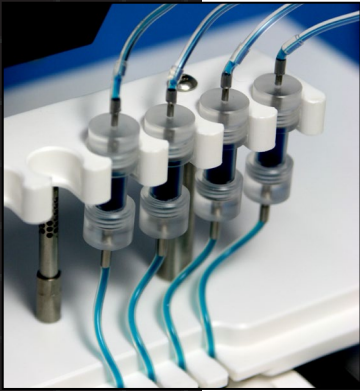

Perifusion Chamber setup

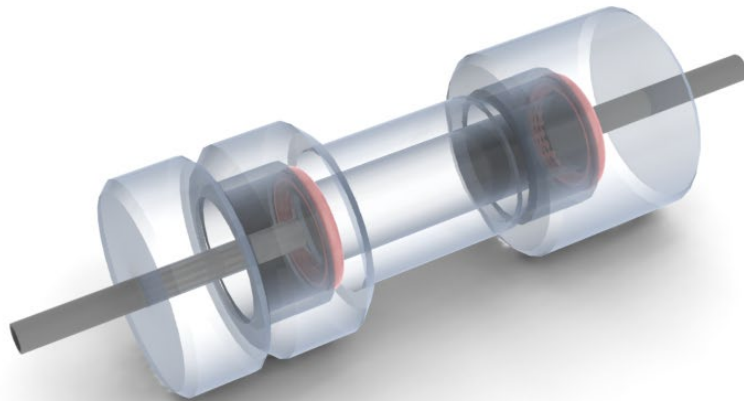


The Perifusion Chamber is specifically designed to use with the Biorep Perifusion System.

It is where all cell reactions will take place. Substances will be delivered through the inlet tube, cells contained in the chamber will react to the substances, and secretions will come out the outlet tube for collection.

- Internal volume = 275 μ L
- Optically clear body, PMMA material
- Easy cell recovery using a micro-centrifuge tube

*Custom sizes and materials available upon request

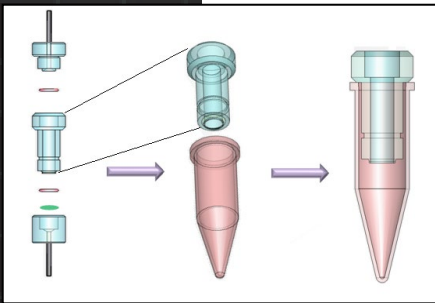

Easy cell recovery